

Haemodialysis

Citrosteril

For Heat Disinfection of Haemodialysis Machines
with Recirculation



Fresenius Medical Care

Citrosteril

For Heat Disinfection of Haemodialysis Machines with Recirculation

Citrosteril – for thermochemical disinfection in haemodialysis machines, e.g. haemodialysis system 4008 or therapy system 5008.

- pH value 1.7 to 2.0
- dissolution of blood residues
- excellent removal of CaCO_3
- disinfection and decalcification in one process
- active ingredients composed of natural substances
- biodegradable
- odourless
- free from colouring additives

Action

The synergistic effect of its components makes Citrosteril a potent disinfectant solution.

Citrosteril at 84°C has a broad spectrum of microbicidal activity and works bactericidal, virus inactivating (HBV, HCV, HIV) and fungicidal.



Specifications and order information:

100g Citrosteril contains:

21g citric acid 1-hydrate; lactic acid, malic acid


Unit	Language combination	Art. No.
1 × 5 L	multilingual	508 533 1
6 × 2 L	multilingual	508 536 1

Literature:

Solbach W, Universität zu Lübeck: Verification of the sporicidal efficacy of the product Citrosteril at 85°C in the Fresenius 5008 dialysis machine, 20.12.2002

Labor Dr. Merk&Kollegen, Ochsenhausen: Antiviral efficacy of Citrosteril against bovine Parvovirus, 9.9.2005

Further information is available on request.

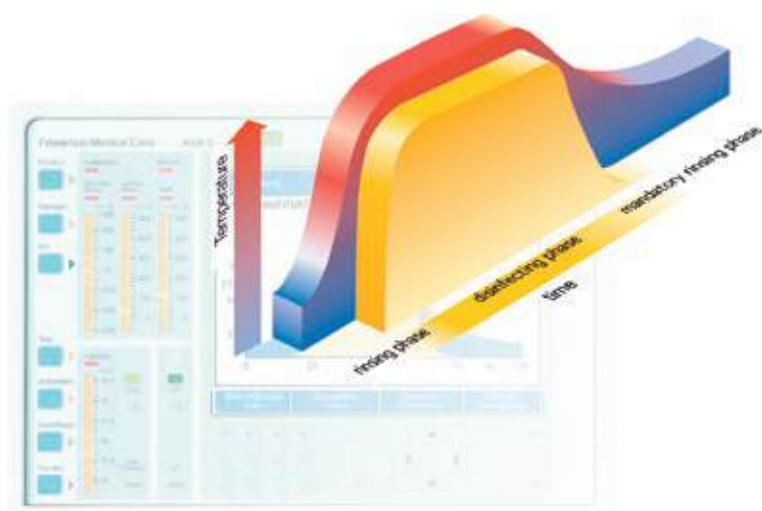

Fresenius Medical Care

Fresenius Medical Care Deutschland GmbH · 61346 Bad Homburg v. d. H. · Germany · Phone: +49 (0) 6172-609-0 · Fax: +49 (0) 6172-609-2191

Head office: Elise-Kröner-Straße 1 · 61352 Bad Homburg v. d. H.

www.fmc-ag.com